

	Pro Skid-E	4BELL STAIR	4BELL CLASS			
Product code	SK10001E	ST10437B	ST10438B			
Height opened	1580 mm	1360 mm	905 mm			
Height folded	1060 mm	930 mm	930 mm			
Length with handles closed	1110 mm	750 mm	750 mm			
Length with handles opened	1340 mm	1170 mm	1170 mm			
Width	550 mm	510 mm	440 mm			
Thickness folded	355 mm	310 mm	310 mm			
Weight	15 kg	14.5 kg	9.3 kg			
Materials	Al, PVC, PE, PU	Steel, Al, PU, PVC	Steel, Al, PVC			
Front wheels	Ø 100 mm	Ø 100 mm	Ø 100 mm			
Rear wheels	Ø 200 mm	Ø 150 mm	Ø 150 mm			
Sheet colours	Black	Yellow, black	Yellow, black			
Load capacity	150 kg	180 kg	180 kg			
Belts polypropylene 50 mm straps	2	2	2			
Welded frame	YES	NO	YES			
SPEN-TEX sheet	YES	YES	YES			
Seat padding	YES	YES	YES			
Other serie codes	Skid-Ok Skid-Ok / B with armrests Skid-E Skid-E Ready fixed backrest Pro Skid-E silver	SK20003E SK20002E SK20001E SK20101C SK10201E	4BELL Air 4BELL Comfort 4BELL Class	ST10439B ST10436B ST10438B	4BELL Air 4BELL Comfort 4BELL Stair	ST10439B ST10436B ST10437B

Are you aware of all **laws** applicable to Spencer chairs?

Spencer chairs series has been designed in compliance with all relevant norms and regulations

	Skid Series	4Bell Series	Compact Series
UNI EN 1865-4		•	
UNI EN 1789 if used with the dedicated fastener		•	•

10 G FASTENERS FOR CHAIRS

Designed for the best safety


Spencer offers a complete range of fasteners for evacuation and transport chairs. Each model is specifically engineered to provide the most secure support to the chair inside an emergency vehicle. The combination of each transport and evacuation chair with its own specific fixing device creates a system able to offer maximum compactness for storage, intuitive use and ideal safety.



4BELL STAIR MAX
10 G fastener for 4BELL STAIR chair
ST10441B



4BELL MAX
10 G fastener for chair 4BELL
ST10432B

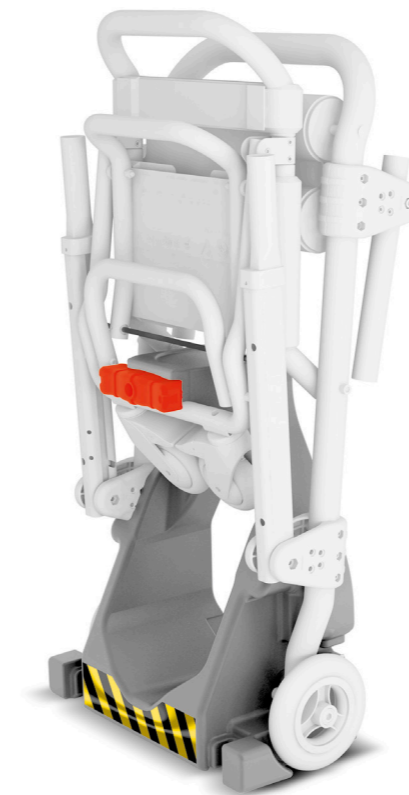

WALL BRACKET
Wall Bracket For Skid/Pro Skid
SK14000B

WALL BRACKET
Wall fastener for SKID-OK
and SKID-OK B evacuation chairs
SK14001B

The value of a choice is measured in the future.

Skid Series / 4BELL Series
Introducing the ultimate expression of Spencer's efficiency, design and innovation.

spencer.it



SKID SERIES

Evacuation stair chair

The simplicity of Spencer's projects is the result of the intuitive nature of its interface. In pursuit of the sublime balance where nothing is added nor removed, **Spencer creates objects and instruments that can be used intuitively**, "without thinking", even in one of the most difficult sectors for designers, that of items for use in evacuation scenarios as well on a daily basis. An absolutely innovative project, part of the exclusive heritage of a unique worldwide manufacture.

Unique features

- PVC sheet easily cleanable equipped with quick release buckles to adjust the tension
- Ergonomic footrest
- Pivoting front wheels with brake
- Back wheels of large diameter to improve the sliding
- Mobile headrest sheet to adapt to the patient's height

Patented integrated system for self-positioning of the front wheels. It allows an automatic correct positioning of the front wheels when the chair is closed for storage.

- SELF-BRAKING SLIDES
- INTUITIVE USE
- MAXIMUM COMFORT



4BELL SERIES

10 G ultralight chairs

Spencer 4BELL chair represents evolution in the concept of evacuation. **Spencer 4BELL is 4 times revolutionary for the world of emergency**. Spencer 4BELL has no screws and the structure is made solid thanks to an innovative system of integrated self-locking elements. This new technique avoids any protrusions dangerous for both the patient and operators and eliminates any risk of instability due to fasteners that tend to loosen or unscrew.

Unique features

- Easyblock system for intelligent opening of front telescopic handles
- PVC sheet easy to disinfect
- Two rear foldable handles to save space and to facilitate the lifting in every situation
- Spencer 4BELL is foldable to save space when in storage
- Self-positioning patented front wheels with pivot and roll brake avoids any hindrance and noise during transportation
- Semirigid seat in SPEN-TEX and closed foam able to offer excellent comfort even over long period of time

Certified 10 G when used together with the dedicated fixing system.

10 G

- SECURE AND RELIABLE
- LIGHTWEIGHT
- ROBUST FRAME





The complete range of 4BELL Evacuation chairs is certified 10 G when combined with the dedicated fixing system 4BELL Max. All models are equipped with two integrated safety belts; one positioned on the backrest, the second on the seat.

4BELL CLASS THE ONLY ONE YOU CAN LIFT WITH A FINGER

4BELL Class weighs only 9.4 kg and it's characterized by strong foldable frame. Secure, compact, and intuitive with Easyblock system, you can open front handles with a touch. 4BELL Class is also equipped with Self-positioning wheels with brakes to facilitate all procedures.

4BELL COMFORT REVOLUTIONARY, ALSO IN COMFORT

4BELL Comfort includes accessories, such as armrests and headrest, to improve patient comfort and stability. Such upgrades encourage correct posture and limits the patient's movements. These accessories also assist the rescuer in his job.

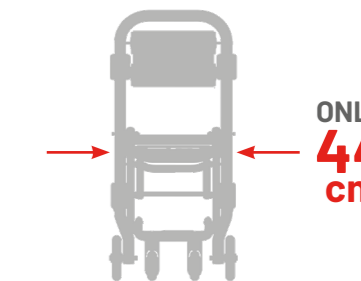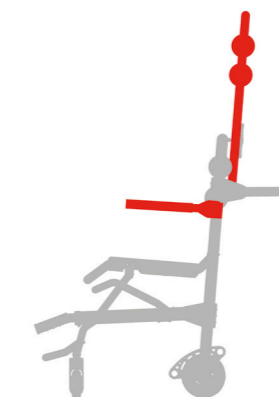

4BELL AIR THE EVACUATION CHAIR THAT'S READY TO FLY

4BELL Air is the narrow version of 4BELL series that is particularly suitable on aircraft, ships, public buildings, shopping malls and old buildings, where standard devices would be too wide.

4BELL STAIR THE EFFORTLESS STAIR CHAIR

4BELL Stair is equipped with slides for stair descent. The headrest design offers to the patient comfort and stability and it is used as a handle during descent. Lightweight, compact, rapid and agile, now descent is an easy task.

