


CARRERA TF PRO

Self-Loading stretcher

With housing for spine board

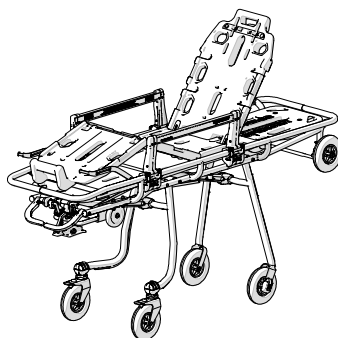

Carrera Pro is a self-loading stretcher for patient transport usable in rescue operations on streets and rescue vehicles. The legs fold away independently, activating the dedicated commands and they open automatically during the unloading procedure. The spinal board can be placed in a dedicated housing located under the patient board.

Specific features

Polyethylene board easy to clean
Color coded commands
Adjustable backrest
Antishock positions Trendelenburg and Fowler
2 fixed wheels and 2 pivoting with large diameter to make movement easier
Brakes on rear wheels
Foldable side bars <i>Safe bar</i>
Spinal board easy to remove, making immobilization procedures faster

Accessories

ST42100A	FP Posterior fastener for self loading stretchers
ST42200A	FA Anterior fastener for self loading stretchers
ST42712A	HUX INTEGRAL FASTENER
EN90003C	END-T Table for stretchers 10g cert.
ST00491A	STX 90 Telescopic headrest for stretchers
ST00497B	DNA STRAP Belt w/integrated re-winding system
ST00498B	DNA STRAP Thorax. Belt w/integr. re-wind. syst.
IF01047C	Track 4-30 telescopic IV holder
IF01049B	Track 5 telescopic IV holder for Carrera
CB09025C	Tanker oxygen tank holder for stretcher side bar
CB09028C	Tanker support for Safe Bar
ST02061B	B-BAK PIN YELLOW SPINE BOARD W/PINS
ST70018A	QMX776 Anatomic mattress

Note Is possible to use only the spinal board Spencer B-Bak PIN not included

Standard Equipment

ST70002A	STX 702 two pieces black reflex metal belts
----------	---

Technical data

Length	1980 ± 10 mm
Width	585 ± 10 mm
Weight	36 ± 1 kg
Loading capacity	170 kg
Wheels	ø200 mm
Height loading carriage wheel	60 cm
Frame material	Stainless steel
Side bars length	680 mm
Side bars Height	200 mm

Class I Med.Device according to Dir. 93/42/CEE

CARRERA TF PRO ALTA T/CROMO

P/NERO

CA01251E



CND Classification

V0899

Registration number

148438

Rev.0 (09/03/2018)

UNCHECKED COPY – further revisions will be available on <http://support.spencer.it>
 Spencer Italia S.r.l. Sala Baganza (PR) Italia Tel. +39.0521.541111 Fax +39.0521.541222

spencer.it