


CINCO MAS

10g Self loading stretcher with 5 heights

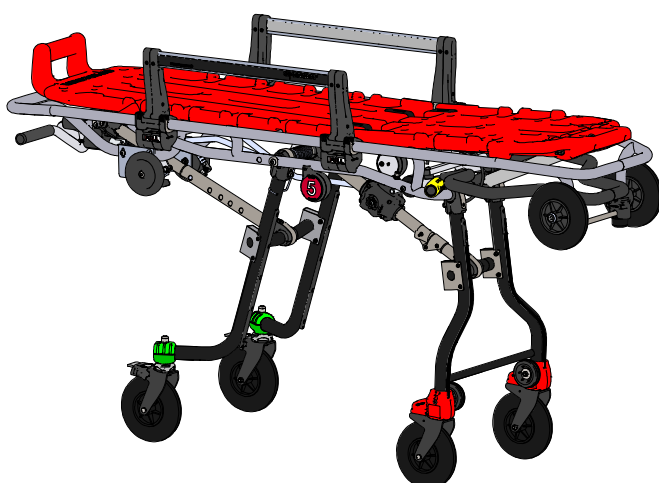

Spencer Cinco Mas is a self-loading stretcher with 5 adjustable heights for patient transport. It is a compact and light stretcher, usable in rescue operations on streets and rescue vehicles. The legs fold away independently, activating the dedicated commands during the unloading procedure, they open automatically. 10G certified if used with the fastener *Sensor Lock*.

Specific features

- Polyethylene board easy to clean
- Adjustable backrest
- Antishock position Trendelenburg
- Extreme maneuverability thanks to the Twist system which, when operated, makes pivoting all 4 wheels
- Brakes on rear wheels
- Foldable side bars *Safe bar*
- 5 Adjustable heights system
- Integrated SHARP SENSOR LOCK system on the leg unlocking handle, which allows to unlock the stretcher from the fastener without acting on the fastener
- Adjustable heights unlocking system logically placed to improve usability
- The patented front leg opening and closing system, reduces the effort required by the operator, and allows to make such movement much more smooth

Needed optional accessories

- CR90013B Loading trolley extender for Cinco Mas

Optional accessories

- ST42706 10G fastener sensor lock (Mandatory for 10g versions)
- EN90003C END-T Table for stretchers 10g cert.
- ST00491A STX 90 Telescopic headrest for stretchers
- ST00499B STX 499 – 4 pcs adjustable thorax belt
- IF01047C Track 4-30 telescopic IV holder
- IF01049B Track 5 telescopic IV holder for Carrera
- CB09025C Tanker oxygen tank holder for stretcher side bar
- CB09028C Tanker support for Safe Bar

Standard Equipment

- | | |
|----------|---|
| ST70000A | QMX 777 BLACK ANATOMIC MATTRESS, |
| ST00497B | DNA STRAP Belt w/integrated re-winding system |
| ST00498B | DNA STRAP Th.Belt w/integr.re-winding system |

Technical data

Length	1970 mm
Width	570 mm
Weight	42 Kg
Loading capacity	250 kg
Wheels	ø200mm ± 5%
Height loading carriage wheel	65 cm
Frame material	Stainless steel
Side bars length	680 mm
Side bars Height	200 mm

Class I Med.Device according to Dir. 93/42/CEE

CINCO MAS 65 T RED BOARD 10G
CERTIFIED

CND Classification
V0899



Rev.0 (12/03/2019)

UNCHECKED COPY – further revisions will be available on <http://support.spencer.it>
Spencer Italia S.r.l. Sala Baganza (PR) Italia Tel. +39.0521.541111 Fax +39.0521.541222

spencer.it