


CINCO MAS

10g Self loading stretcher with 5 heights

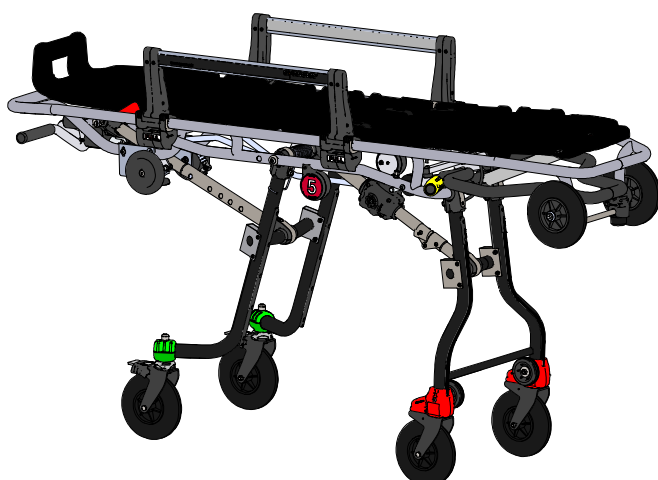

Spencer Cinco Mas is a self-loading stretcher with 5 adjustable heights for patient transport. It is a compact and light stretcher, usable in rescue operations on streets and rescue vehicles. The legs fold away independently, activating the dedicated commands during the unloading procedure, they open automatically. 10G certified if used with the fastener *Sensor Lock*.

Specific features

- Polyethylene board easy to clean
- Adjustable backrest
- Antishock position Trendelenburg
- Extreme maneuverability thanks to the Twist system which, when operated, makes pivoting all 4 wheels
- Brakes on rear wheels
- Foldable side bars *Safe bar*
- 5 Adjustable heights system
- Integrated SHARP SENSOR LOCK system on the leg unlocking handle, which allows to unlock the stretcher from the fastener without acting on the fastener
- Adjustable heights unlocking system logically placed to improve usability
- The patented front leg opening and closing system, reduces the effort required by the operator, and allows to make such movement much more smooth

Needed optional accessories

- | | |
|----------|--|
| CR90013B | Loading trolley extender for Cinco Mas |
|----------|--|

Optional accessories

- | | |
|----------|--|
| ST42706C | 10G fastener sensor lock ^(Mandatory) |
| EN90003C | END-T Table for stretchers 10g cert. |
| ST00491A | STX 90 Telescopic headrest for stretchers |
| ST00499B | STX 499 – 4 pcs adjustable thorax belt |
| IF01047C | Track 4-30 telescopic IV holder |
| IF01049B | Track 5 telescopic IV holder for Carrera |
| CB09025C | Tanker oxygen tank holder for stretcher side bar |
| CB09028C | Tanker support for Safe Bar |

Standard Equipment

- | | |
|----------|---|
| ST70000A | QMX 777 BLACK ANATOMIC MATTRESS, |
| ST00497B | DNA STRAP Belt w/integrated re-winding system |
| ST00498B | DNA STRAP Th.Belt w/integr.re-winding system |

Technical data

Length	1970 mm
Width	570 mm
Weight	43 Kg
Loading capacity	250 kg
Wheels	ø200mm ± 5%
Height loading carriage wheel	75 cm
Frame material	Stainless steel
Side bars length	680 mm
Side bars Height	200 mm

Class I Med.Device according to Dir. 93/42/CEE

CINCO MAS 75 T BLACK BOARD
10G CERTIFIED

CND Classification
V0899



Rev.0 (12/03/2019)

UNCHECKED COPY – further revisions will be available on <http://support.spencer.it>
Spencer Italia S.r.l. Sala Baganza (PR) Italia Tel. +39.0521.541111 Fax +39.0521.541222

spencer.it