


## POP

### Oxygen distribution line AFNOR

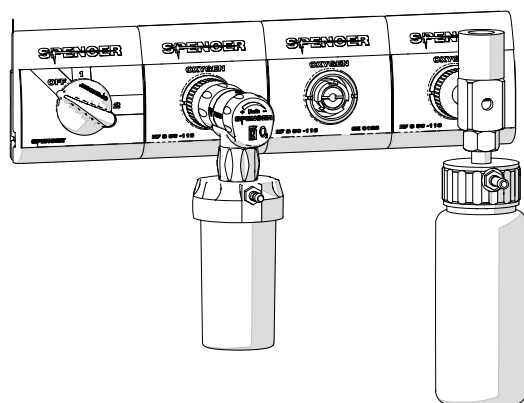

System for wall recessed mounting, with three medical oxygen output units, to be installed inside a medical vehicle on dedicated surfaces.

#### Specific features

- Simple and safe use guaranteed by the adoption of gas-specific connections and quick release mechanism
- Smooth surfaces for an easy cleaning
- The locking mechanism, allows easy installation
- Equipped with flowmeter complete of humidifier and with venturi suction unit
- Complete of 2 pressure reducers and manual tank exchanger

#### Technical data

Overall dimensions (without connected devices)	400x100x86 ±5 mm
Dimensions of the non-recessed part	400x100x36,5 ±5 mm
Connection standard	AFNOR
Materials	Nylon, Al, Brass
Flowmeter steps	From 0 to 15 l/min
Humidifier volume	Max 130ml
Venturi consumption	Max 47 l/mn
Venturi max. flow	25 l/min
Venturi collection canister	300 ml
Operating pressure	3,5-5 bar

#### Compliant with

UNI EN ISO 9170  
NF S 90-116  
EN 1789

Ass. Med. Dev. According to Dir. 93/42/CEE  
MOD.G012 DISTRIBUTION LINE  
FOR MEDICAL GAS  
OX04392C



Rev.1 (15/01/2018)

UNCHECKED COPY – further revisions will be available on <http://support.spencer.it>  
Spencer Italia S.r.l. Sala Baganza (PR) Italia Tel. +39.0521.541111 Fax +39.0521.541222  
[www.spencer.it](http://www.spencer.it)