


WOW

Specific features

Carrying sheet with 10 handles

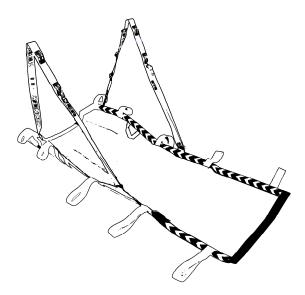

WOW injury transportation sheets are devices to be used to move and transport patients in a lying or seated position.

openio idatai di	
Manufactured with fireproof materials	
Easy to clean thanks to the use of non-absorbent fabrics	
Reinforced nylon webbing to ensure high resistance and durability	
The width of the tape used for the handles, ensures are easy grabbing increasing safety	
On the patient side there are handles dedicated to the patient that make him more collaborative and increase safety	
The special shape, offers high lateral containment during transport	
The supplied adjustable shoulder straps ensure a better distribution of the load on the rescuer	
The longitudinal aluminium rods are removable through	

	limit the flexion of the device.				
	The device also allows the transport of the patient in				
	sitting position				
Requires low storage space					

the dedicated holes and can be folded into 3 parts. When inserted, they increase the rigidity of the sheet, allowing

to distribute the load over a wider surface and allowing to

Reference standards	
	UNI EN 1865-1

Technical data		
Length	1870 ± 10 mm	
Upper part width	960 ± 5 mm	
Lower part width	700 ± 5 mm	
Dimensions when stored (1)	500 x 600 x 150 mm	
Dimensions of single folded rod	74 x 450 x 30 ± 5 mm	
# Handles	10	
Strap width	50 mm	
Length of upper handles	160 mm	
Length of lower handles	260 mm	
Rod length	1210 ± 10 mm	
Number of side strap	6	
adjustments	O	
Capacity	150 kg	
Weight	4,2 ± 0,2 kg	

(1) The measurement refers to wrapped sheets. This is purely indicative, as the device can be folded or wrapped in numerous ways

Class I MD compliant with UE Reg. 2017/	745
WOW 10 HANDLES TRANSFER	
SHEET	((
QM40601B	
CND Classification	V0899