


# MSF

## Field stretcher

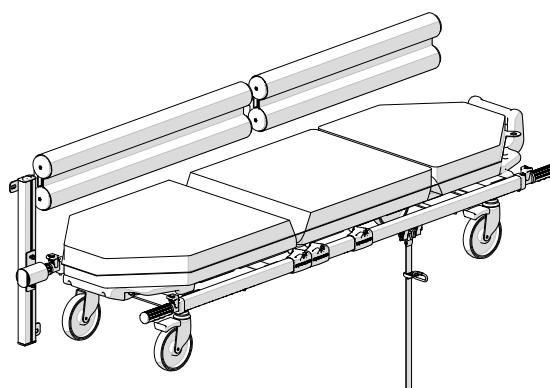

PRIME series stretchers MSF are devices for lifting and carrying people for use in field rescue situations. Parking is permitted preferably when using the parking stops installed in designated locations within the rescue area.

### Specific features

- The installation can be done using the original fixing elements and fasteners of the vehicle, without requiring new holes or additional processing
- The system locks both the front and rear part of the stretcher. The stretcher is released from the vehicle opening the rear door.
- Backrest for patient who do not require to be transported in lying position or for rescuers
- Polyethylene patient board easy to clean
- Backrest adjustable in 4 positions
- Front and rear telescopic handles
- Wheels have large diameter to increase sliding action
- Brakes on outer wheels and quick locking belt to stabilize the device when the vehicle is in motion
- Soft mattress with 3 segments to increase patient comfort
- Included 2 belts with plastic hook for restrain the patient

### Technical data

Stretcher length	1895
Stretcher width	575
Possible anchorage points	3
Materials	Al, PVC, PE, Nylon, Steel
Overall weight	25 ± 2kg

### Accessories

IF01043B	IV holder Track 4
ST00497B	DNA Strap with integrated re-winding system
ST00498B	DNA Strap thoracic belt with re-winding system
ST00499B	STX499 – Thorax belts

Class I MD compliant with UE Reg. 2017/745

MOD.026 KIT FOR STRETCHER FOR  
TOYOTA HZJ78

ST00052

CND Classification



V0899

Rev.0 (03/08/2021)

UNCHECKED COPY – further revisions will be available on <http://support.spencer.it>  
Spencer Italia S.r.l. Sala Baganza (PR) Italia Tel. +39.0521.541111 Fax +39.0521.541222

[www.spencer.it](http://www.spencer.it)